

YOUR WEIGHT SCALE



1 - Display 2 - Measurement Switch

SCALE DISPLAY



GETTING STARTED

1

Place the scale on a solid level surface, such as tile, wood, or laminate flooring. If the scale must be used on a carpet, be sure to attach the plastic add-on feet.

2



Open the battery cover. Insert “AA” batteries and replace the battery cover.

3



Set the slide switch to the “lb” setting to ensure measurement is in pounds.


4



Press the measurement switch gently.

5




All display segments are displayed for several seconds. Wait until the  symbol and “0.0 kg /0.0 lb” are displayed.

6



Step on the scale gently and stand still during the measurement.

7

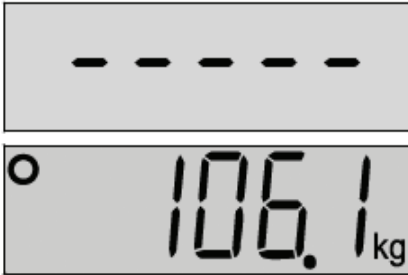

Your weight is displayed after the symbol is displayed. After about fifteen seconds, the scale will automatically turn off.